



Contraindications:

Absolute (When the product should not be used at all)

- Demonstrated sensitivity to hydrocolloid (rare)

Relative (When the product can be used after healing)

- Unhealed wound at application site (wait to heal)
- Active rash at application site (wait to heal)
- Active inflammation or infection

For use on intact skin.

If redness or irritation occurs, stop use and consult a healthcare practitioner.

Any serious incident while using this device should be reported to the manufacturer and competent authority.



Manufactured by: **BioDerm, Inc.**
12320 73rd Ct. N., Largo, FL 33773, USA
Orders: (800) 814-3174 Support: (800) 864-6755
CustomerCare@BioDermInc.com
www.BioDermInc.com



NOT MADE WITH NATURAL RUBBER LATEX

**Not included in part numbers ending with NA*

©BioDerm, Inc. All rights reserved. BioDerm® and CathGrip® are registered trademarks of BioDerm, Inc. International European patent EP(UK/NL) 2501429.

Application Tips



Make sure the patient's skin is **COMPLETELY DRY** before applying CathGrip®.



Cleansing wipes and Chlorhexidine can cause excess moisture underneath the skin that can reduce CathGrip wear time. Make certain that the application site for CathGrip is as dry and moisture free as possible, before applying.



Wear time. The first application for each patient will last shorter than the subsequent applications of the product. This is due to the hydrocolloid material that absorbs underlying moisture in the skin. CathGrip can be worn for up to 24 hours.



Single-Use. Once removed, dispose of in a trash can. Product is not reusable and could cause product malfunction or risk of infection if reused.

Instructions

For Best Results Follow These Steps

• Prep • Peel • Press • Place • Pull • Remove

Contents: • Hydrocolloid Universal Securement Device

1 Preparation


Wipe & Dry

Clean entire application site with a non-moisturizing soap. Allow to dry.

Hydrocolloid
Adhesive Oval

Soft Strap
with No-Slip
Grip

Release
Papers

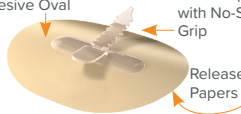


2 Peel




Remove the release papers from the hydrocolloid adhesive oval.

3 Press

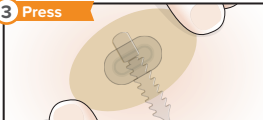


Press the hydrocolloid oval to the skin for 15 seconds to seal.

4 Place

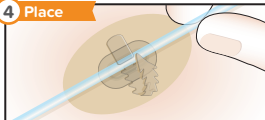


Place the tube on the soft, no-slip grips.

5 Pull & Secure

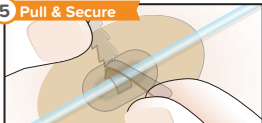

Pull the straps through to secure.

6 Remove

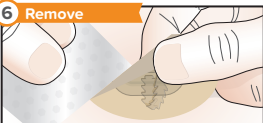

Remove with soapy water or per facility's removal procedures.

REMOVAL: Remove and replace CathGrip when the adhesive turns a white color, when it is loose, or when it is time to remove.

ATTENTION: Following application, leave a one inch space around the CathGrip® device when using any cleansing products, creams or lotions. Contact with the hydrocolloid will result in early dislodgement of the securement device.